

# **Extremity Whirlpools**

Division of Acorn Engineering Company









# E Series 22 Gallon Stationary Whirlpool with Pedestal

**E-22-SP** 



E-22-SP

Whitehall E-22-S Stationary Whirlpool with Pedestal is recommended to facilitate hydrotherapy of the foot, ankle and leg, to

hydrotherapy of the foot, ankle and leg, to below the knee, with better coverage and greater patient comfort.

Tank Capacity is 22 gallons.

Whirlpool Tank is 28"L x 15"W x 18"D and is seamless welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. Tank rim is reinforced with a concealed stainless steel rod. Tank bottom is seamless, coved (rounded) design for effective cleaning. Tank is fitted with a fixed raised pedestal for mounting to floor.

A 2" Combination Drain Overflow and a Filler Spout are provided.

**Turbine Raising and Lowering Device** permits adjustment of desired height and direction of water agitation towards the area of the body. It is spring balanced and provided with a locking device.

**Turbine Assembly** is equipped with Underwriters Laboratory listed 1/2 HP jet pump motor protected with an ABS plastic shell, and automatic thermal overload protector and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze.

**Thermometer** is provided with a dual scale, adjustable calibration, and a stem inserts into a retainer to prevent rattle during whirlpool agitation.

### **GUIDE SPECIFICATIONS**

Provide Whitehall Stationary Whirlpool (specify model number and options) with a 22-gallon capacity, 28"L x 15"W x 18"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. Tank shall have a reinforced rim and shall be fitted with a fixed raised pedestal. Bottom of tank shall be a seamless, rounded coved design. Provide a 2" combination Drain Overflow, and a filler spout. Provide a turbine assembly with raising and lowering device which functions as an agitator and permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP jet pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a thermometer with a stem retainer to prevent rattle when whirlpool is in operation.

Please visit www.whitehallmfg.com for most current specifications.





# **Extremity Whirlpools**

### **MODEL NUMBER AND OPTION SELECTIONS**

### Series

☑ E Extremity Series

### **Capacity**

✓ -22 Gallons

### **Type**

✓ -SP Stationary Whirlpool with Pedestal

### **Turbine Electrical Requirements**

- ☐ 115VAC / 50Hz
- ☐ 115VAC / 60HZ
- ☐ 230VAC / 50Hz
- ☐ 230VAC / 60HZ

Country of Use

### For Mixing Valve:

Refer to supplemental MXT15-OTG and DU-1 sheet for details.

**Suffix Options** (See Whirlpool Options Section for Details)

Enviro-Glaze® Color Powder Coating -□ -EG

Specify Color:

https://www.whitehallmfg.com/enviro-glaze

🔲 -GFI Ground Fault Interrupter

☐ -LOGO Applied Logos: Team, Facility, Company

Separate Drain and Overflow Assembly ☐ -SDO

## **ELECTRIC REQUIREMENTS:**

Amps: 6.9/6.2 Hertz: 50/60 Volts: 115 GFI Receptacle


### **ELECTRIC REQUIREMENTS:** Amps: 3.5/3.1 Hertz: 50/60 Volts: 230 GFI Receptacle


Certified per US & Canadian Medical Standards: UL 60601-1 & CSA C22.2 No

Certified per European Standards; **(**EC 60601-1





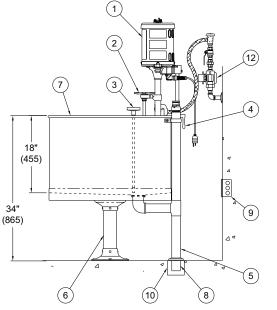

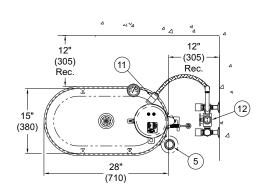

### NOTES:

- 1. 1/2" HP Jet Pump Motor
- 2. Variable Pressure Control
- 3. Thermometer
- 4. Turbine Raising & Lowering Device

### STATIONARY E-22-SP SERIES WHIRLPOOL

- 5. Combination Drain & Overflow Assembly
- 6. Mounting Pedestal
- 7. Stainless Steel Tank
- 8. 2" O.D. x 5" Tailpiece
- 9. 3 Pole GFI Receptacle by Others
- 10. Floor Drain by Others
- 11. Filler Spout
- 12. Mixing Valve, Shown For Reference Only

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility fo use of void or superseded data © Whitehall Manufacturing, Member of Morris Group International.

Please visit www.whitehallmfg.com for most current specifications

▲ WARNING: Cancer and Reproductive Harm - www.P65Warnings.ca.gov

SELECTION SUMMARY & APPROVAL FOR MANUFACTURING		
Model Number & Options		Quantity
Company		Date
Contact	Title	
Approval for Manufacturing/Signature		

Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 591-9280 • (626) 968-6681