

Extremity Whirlpools

Division of Engineering Company®



Jse



E Series 27 Gallon Mobile Whirlpool

E-27-M



Whitehall E-27-M Mobile Whirlpools are recommended to facilitate hydrtherapy of the upper extremities and the lower leg below the knee.

Tank Capacity is 27 gallons.

Whirlpool Tank is 28"L x 15"W x 21"D and is seamless welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. Tank rim is reinforced with a concealed stainless steel rod. Tank bottom is a seamless, coved (rounded) design for effective cleaing.

Heavy Duty Swivel Casters support all mobile models. Rear pair of casters are lockable for stabilizing during Whirlpool use.

Turbine Raising and Lowering Device permits adjustment of desired height and direction of water agitation towards the area of the body during Whirlpool use. It is spring balanced and provided with a locking device.

Turbine Assembly functions both as agitator and emptying device. It is equipped with Underwriters Laboratory listed 1/2 HP jet pump motor protected with an ABS plastic shell, automatic thermal overload protector, and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze.

Drain option -SDP has a bottom drain, leading to a separate, concealed water evacuation pump. A hose is attached to the side of the unit to facilitate the evacuation.

Thermometer is provided with a dual scale, adjustable calibration, and a stem inserts into a retainer to prevent rattle during whirlpool agitation.

GUIDE SPECIFICATIONS

Provide Whitehall Mobile Whirlpool (specify model number and options) with a 27-gallon capacity, 28"L x 15"W x 21"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. Bottom of tank shall be a seamless. rounded coved design. Whirlpool shall have heavy duty swivel casters with locking device on rear pair of casters. Provide a turbine assembly with raising and lowering device which functions both as agitator and emptying device, permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP jet pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a thermometer with a stem retainer to prevent rattling when whirlpool is in operation.

Please visit www.whitehallmfg.com for most current specifications.





Extremity Whirlpools

MODEL NUMBER AND OPTION SELECTIONS

Series

☑ E Extremity Series

Capacity

✓ -27 Gallons

Type

₫ -M Mobile Whirlpool

Turbine Electrical Requirements

☐ 115VAC / 50Hz ☐ 115VAC / 60HZ ☐ 230VAC / 50Hz ■ 230VAC / 60HZ

Country of Use

For Mixing Valve:

Refer to supplemental MXT15-OTG and DU-1 sheet for details.

Suffix Options (See Whirlpool Options Section for Details)

Auxiliary Drain Hose Assembly 🔲 -ADH

Enviro-Glaze[®] Color Powder Coating -□ -EG


Specify Color:

https://www.whitehallmfg.com/enviro-glaze

🔲 -GFI Ground Fault Interrupter

-LOGO Applied Logos: Team, Facility, Company

☐ -SDP Separate Drain Pump

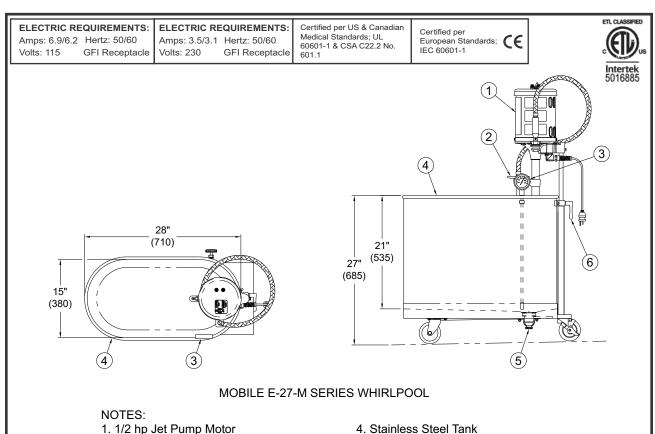

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility fo use of void or superseded data © Whitehall Manufacturing, Member of Morris Group International.

2. Variable Pressure Control

3. Thermometer

Please visit www.whitehallmfg.com for most current specifications

5. Auxiliary Drain Valve

6. Turbine Raising & Lowering Device

MARNING: Cancer and Reproductive Harm - www.P65Warnings.ca.gov

<u> </u>	'	0	0
SELECTION SUMMARY & APPROVAL FOR MANUFACTURING			
Model Number & Options	Quantity		
Company	Date		
Contact	_Title		
Approval for Manufacturing/Signature			

Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 591-9280 • (626) 968-6681