

## **Lo-Boy Whirlpools**

Division of Acorn Engineering Company









# L Series 105 Gallon Stationary Whirlpool

L-105-S



Whitehall L-105-S Stationary Whirlpool is recommended to facilitate hydrotherapy of extremities, hip and back area and allows spaces for exercising.

Whirlpool Tank is 66"L x 24"W x 18"D and is seamless welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. Tank rim is reinforced with a concealed stainless steel rod. Tank bottom is seamless, coved (rounded) design for effective cleaning. Tank is also provided with an extended stainless steel skirt flush to the floor line to conceal internal drain piping. A 2" Combination Drain Overflow and a Filler Spout are provided.

**Turbine Raising and Lowering Device** permits adjustment of desired height and direction of water agitation towards the area of the body during Whirlpool use. It is spring balanced and provided with a locking device.

**Turbine Assembly** is equipped with Underwriters Laboratory listed 1/2 HP jet pump motor protected with an ABS plastic shell, automatic thermal overload protector, and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze.

**Thermometer** is provided with a dual scale, adjustable calibration, and a stem inserts into a retainer to prevent rattle during whirlpool agitation.

#### **GUIDE SPECIFICATIONS**

Provide Whitehall Stationary Whirlpool (specify model number and options) with a 105-gallon capacity, 66"L x 24"W x 18"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. Tank shall have a reinforced rim and shall be fitted with a stainless steel skirt flush to the floor line. Bottom of tank shall be a seamless, rounded coved design. Provide a 2" combination Drain Overflow and a filler spout. Provide a turbine assembly with raising and lowering device which functions as an agitator and permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP jet pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a thermometer with a stem retainer to prevent rattle when whirlpool is in operation.

Please visit www.whitehallmfg.com for most current specifications.





# **Lo-Boy Whirlpools**

#### **MODEL NUMBER AND OPTION SELECTIONS**

#### Series

✓ L Lo-Boy Series

#### **Capacity**

₫ -105 Gallons

#### Type

Stationary Whirlpool

#### **Turbine Electrical Requirements**

- ☐ 115VAC / 50Hz
- ☐ 115VAC / 60HZ
- ☐ 230VAC / 50Hz ☐ 230V
  - ☐ 230VAC / 60HZ

Country of Use

### For Mixing Valve:

Refer to supplemental MXT15-OTG and DU-1 sheet for details.

**Suffix Options** (See Whirlpool Options Section for Details)

- -AS3 Adjustable Suspension Seat
- ☐ -EG Enviro-Glaze® Color Powder Coating -

Specify Color:\_

https://www.whitehallmfg.com/enviro-glaze

- ☐ -GFI Ground Fault Interrupter
- ☐ -HR Head Rest, Adjustable
- ☐ -LOGO Applied Logos: Team, Facility, Company

Logo

- -RS Reclining Seat
- -SDO Separate Drain Overflow Assembly
- ☐ -TTS2 Tank Top Seat


## ELECTRIC REQUIREMENTS:

Amps: 6.9/6.2 Hertz: 50/60 Volts: 115 GFI Receptacle ELECTRIC REQUIREMENTS:
Amps: 3.5/3.1 Hertz: 50/60
Volts: 230 GFI Receptacle

Certified per US & Canadian Medical Standards; UL 60601-1 & CSA C22.2 No.

Certified per European Standards; **CE** IEC 60601-1



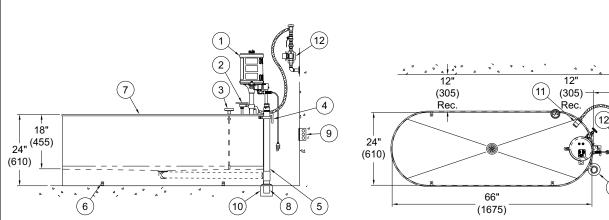

#### NOTES:

- 1. 1/2 HP Jet Pump Motor
- 2. Variable Pressure Control
- 3. Thermometer
- 4. Turbine Raising & Lowering Device

#### STATIONARY L-105-S SERIES WHIRLPOOL

- Combination Drain & Overflow Assembly
- 6. "L" Brackets (4) Four Each
- 7. Stainless Steel Tank
- 8. 2" O.D. x 5" Tailpiece
- 9. 3 Pole GFI Receptacle by Others
- 10. Floor Drain by Others
- 11. Filler Spout
- Mixing Valve, Shown For Reference Only

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility fo use of void or superseded data © Whitehall Manufacturing, Member of Morris Group International.

Please visit **www.whitehallmfg.com** for most current specifications.

▲ WARNING: Cancer and Reproductive Harm - www.P65Warnings.ca.gov

| SELECTION SUMMARY & APPROVAL FOR MANUFACTURING |          |
|--|----------|
| Model Number & Options                         | Quantity |
| Company  | Date     |
| Contact  | Title    |
| Approval for Manufacturing/Signature           |          |

Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 591-9280 • (626) 968-6681