

### **Podiatry Whirlpools**

Division of Acorn Engineering Company



A MEMBER OF





# P Series 22 Gallon Mobile Whirlpool

P-22-M, P-22-MU



P-22-M

Whitehall P-22-M Mobile Whirlpools are recommended to facilitate hydrotherapy of the foot and ankle area, and when provided with a raised undercarriage, may also be used on the upper extremities.

Tank Capacity is 22 gallons.

**P-22-MU** is a Mobile Whirlpool with Undercarriage of heavy duty construction and fixed height which raises the tank for use on the upper extremities.

**Whirlpool Tank** is 28"L x 15"W x 18"D and is seamless welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. Tank rim is reinforced with a concealed stainless steel rod. Tank bottom is a seamless, coved (rounded) design for effective cleaning.

**Heavy Duty Swivel Casters** support all mobile models. Rear pair of casters are lockable for stabilizing during Whirlpool use.

**Turbine Raising and Lowering Device** permits adjustment of desired height and direction of water agitation towards the area of the body during Whirlpool use. It is spring balanced and provided with a locking device.

**Turbine Assembly** functions both as agitator and emptying device. It is equipped with Underwriters Laboratory listed 1/2 HP jet pump motor protected with an ABS plastic shell, automatic thermal overload protector, and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze.

**Drain option -SDP** has a bottom drain, leading to a separate, concealed water evacuation pump. A hose is attached to the side of the unit to facilitate the evacuation

**Thermometer** is provided with a dual scale, adjustable calibration, and a stem inserts into a retainer to prevent rattle during whirlpool agitation.

#### **GUIDE SPECIFICATIONS**

Provide Whitehall Mobile Whirlpool (specify model number and options) with a 22-gallon capacity, 28"L x 15"W x 18"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. Bottom of tank shall be a seamless, rounded coved design. Whirlpool shall have heavy duty swivel casters with locking device on rear pair of casters. Provide a turbine assembly with raising and lowering device which functions both as agitator and emptying device, permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP jet pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a thermometer with a stem retainer to prevent rattling when whirlpool is in operation.

Please visit www.whitehallmfg.com for most current specifications.





## **Podiatry Whirlpools**

#### **MODEL NUMBER AND OPTION SELECTIONS**

#### Series

✓ P Podiatry Series

#### **Capacity**

**☑** 22 Gallons

#### Type

-M Mobile Whirlpool

☐ -M-MU Mobile Whirlpool with Undercarriage

#### **Turbine Electrical Requirements**

☐ 115VAC / 50Hz

☐ 115VAC / 60HZ

☐ 230VAC / 50Hz □ 230VAC / 60HZ Country of Use

For Mixing Valve:

Refer to supplemental MXT15-OTG and DU-1 sheet

**Suffix Options** (See Whirlpool Options Section for Details)

-ADH Auxiliary Drain Hose Assembly

Enviro-Glaze<sup>®</sup> Color Powder Coating -□ -EG

Specify Color:

https://www.whitehallmfg.com/enviro-glaze

🔲 -GFI Ground Fault Interrupter -HL Hydro Lift (Refer To Page HL)

☐ -LOGO Applied Logos: Team, Facility, Company

☐ -SDP Separate Drain Pump

for details.


ELECTRIC REQUIREMENTS: Amps: 6.9/6.2 Hertz: 50/60 Volts: 115 **GFI** Receptacle


**ELECTRIC REQUIREMENTS:** Amps: 3.5/3.1 Hertz: 50/60 Volts: 230 **GFI** Receptacle

Certified per US & Canadian Medical Standards: UL 60601-1 & CSA C22.2 No.

European Standards; (E



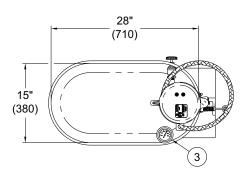

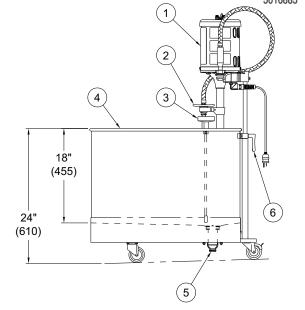

#### MOBILE P-22-M SERIES WHIRLPOOL

### NOTES:

- 1. 1/2 hp Jet Pump Motor
- 2. Variable Pressure Control
- 3. Thermometer

- 4. Stainless Steel Tank
- 5. Auxiliary Drain Valve
- 6. Turbine Raising & Lowering Device

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility fo use of void or superseded data © Whitehall Manufacturing, Member of Morris Group International.

Please visit www.whitehallmfg.com for most current specifications.

▲ WARNING: Cancer and Reproductive Harm - www.P65Warnings.ca.gov

| SELECTION SUMMARY & APPROVAL FOR MANUFACTURING |       |          |  |
|------------------------------------------------|-------|----------|--|
| Model Number & Options                         |       | Quantity |  |
| Company                                        |       | Date     |  |
| Contact                                        | Title |          |  |
| Approval for Manufacturing/Signature           |       |          |  |