

## Surgical Scrub Sink Model Options

Please visit [www.whitehallmfg.com](http://www.whitehallmfg.com) for most current specifications.

**-DTD**

### Digital Timer Display -

This LED display is added to the control panel and displays both minutes and seconds of washing time elapsed, allowing users to monitor the duration of scrub up.




**-DTD**

**-EFW**

### Eye / Face Wash -

A Swing Down Eye and Face Wash unit is provided within the Scrub Sink basin against the back wall to cleanse face and eyes from contaminants. Device is kept in upright position while not in use and activates automatically when swung down. When specified, an ET71-1 Emergency Mixing Valve will be provided and factory pre-piped.

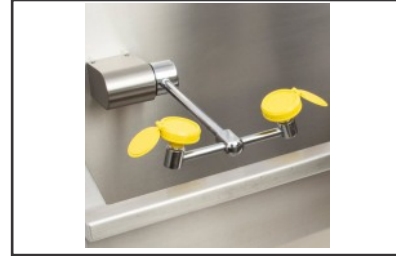


**-EFW**

**-ETC**

### Electronic Timing Controller -

An accurate valve controller provided standard with two scrub times: Short (5 minutes) and Long (10 minutes). It includes a two minute rinse indicator warning light, and a stop cycle button. The scrub cycle times can also be adjusted to fit user's specific protocol.

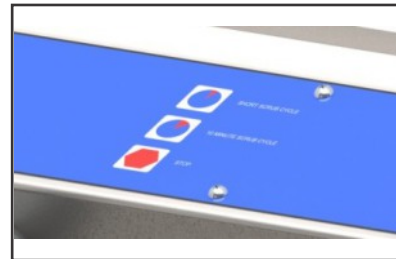

**-ETC**

**-LFN**

### Laminar Flow Nozzle -

Nozzle provides a 1.5 gpm faucet attachment for our Gooseneck spouts. This attachment reduces water flow and provides a smooth, softer, clear, solid stream of water that reduces splashes and helps minimize the spread of airborne bacteria.



**-LFN**

**-PB**

### Pedestal Base -

Pedestal Base provides on-the floor mounting for support and lessens the amount of site or wall prep for mounting and installation. The depth of the Base provides a users with a toe space.



**-PB**